

COVID Lifesaver Mask Manual



Adapter Version
3.0

INDEX



1x

Decathlon Subea Easybreath sizes S/M or M/L



1x

HME filter



1x

Adapter version 3.0
FDM 3D print



The adapter and facemask are meant for personal use. The HME filter needs to be discarded after each uses.

CHECK 3D PRINT






Before assembling the COVID Lifesaver Mask check the FDM 3D printed adapter for loose filament inside or printing errors. Remove any loose filament or, in the case of a printing error, discard the adapter and reprint it.

CONNECT FILTER

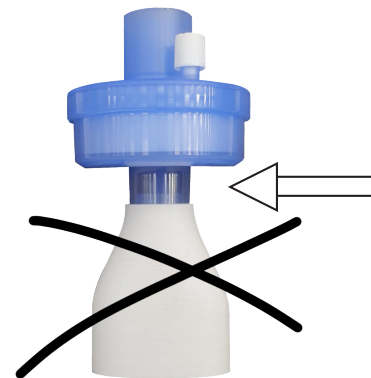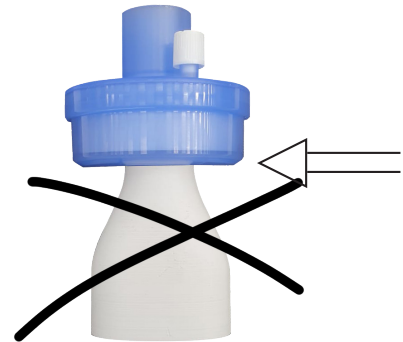

Connect the HME filter to the adapter. Make sure that the filter is pushed in all the way, but be careful not to push too hard; otherwise, the crosspiece inside the filter can break. The filter should stick out about 5mm.

CONNECT ADAPTER

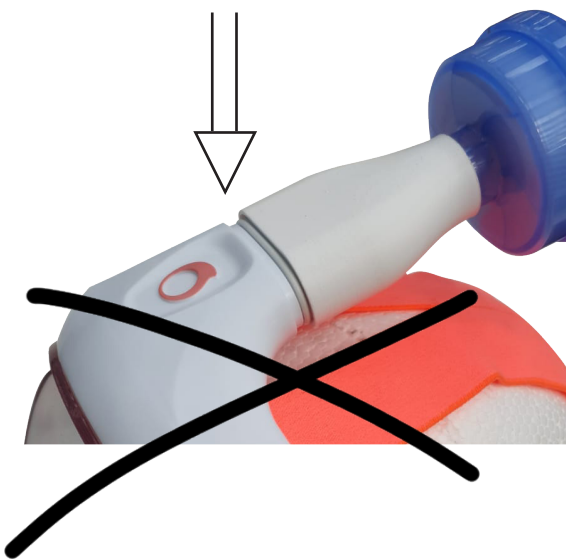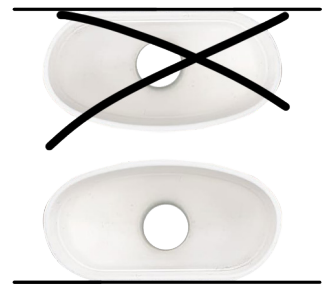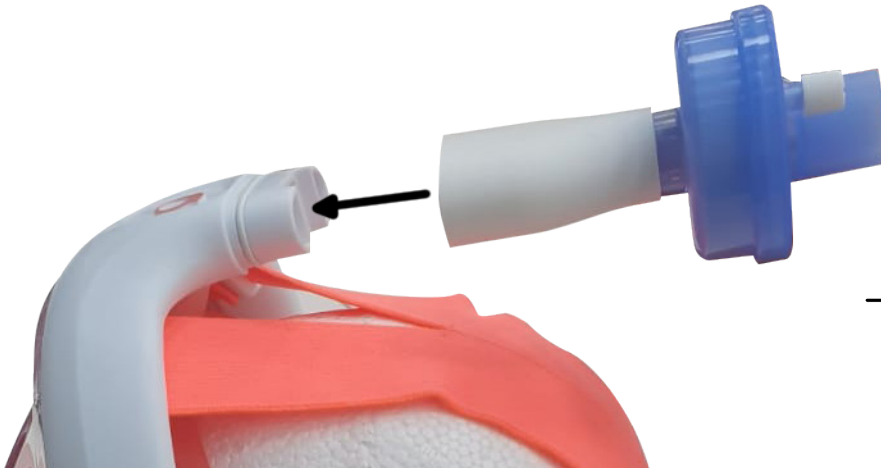

Connect the adapter to the facemask. Make sure the rounded side of the adapter is facing upwards, allowing the adapter to slide on all the way. The adapter should slide on easily without the use of much force.